

SHUR/Trim™

Paraffin Block Dewaxer



TBS®


SHUR/Trim™ . . . Eliminate cut fingers.

An efficient appliance for removing excess paraffin from tissue cassettes following embedding. Simply slide paraffin blocks across the heated grooved surface. Melted paraffin is shaved away to disposable receptacle below. Preparing embedded tissue blocks for insertion into microtome quick-release clamps has never been easier or safer.



Ordering Information:

Cat #	U/M	Description
ST-120	ea	SHUR/Trim™ Paraffin Block Dewaxer; WxDxH (6x9x7)in, (15x23x18)cm; 5lbs, 3kg; 115 VAC, 50/60 Hz, 4 Amps.
ST-PCH	ea	Paraffin Drain Cup Holder.
ST-PDC	pk/cs	Paraffin Drain Cups; 10 ea/pk; 4pk/cs.

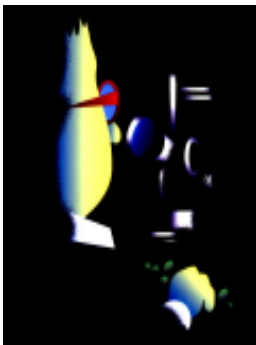

Please visit our website for additional product information, technical bulletins, career opportunities, MSDS files, etc.

www.trianglebiomedical.com

Distributed By:

TBS® TRIANGLE BIOMEDICAL SCIENCES, INC.
3014 CROASDAILE DRIVE • DURHAM, NORTH CAROLINA 27705-2507 • USA • 919.384.9393 • FAX:919.384.9595 • EMAIL: tbs@trianglebiomedical.com

Specifications are current at time of publication; however, they are subject to change without notification. Please contact our North Carolina office or one of our distributors for current prices and availability. SHUR/Trim™ is a trademark and TBS® is a registered trademark of Triangle Biomedical Sciences, Inc., Durham, North Carolina, USA.

Litho in U.S.A. © September 2008

Bulletin # 3.1